


RIESTER RI-FOCUS LED HEADLIGHT

Product Code: UB-RT6090/1

**WIRELESS HIGH
PERFORMANCE LED
HEADLIGHT POWERED BY
LI-BATTERIES (UB-RT6090)**

OR

**RECHARGEABLE
BATTERIES (UB-RT6091)**

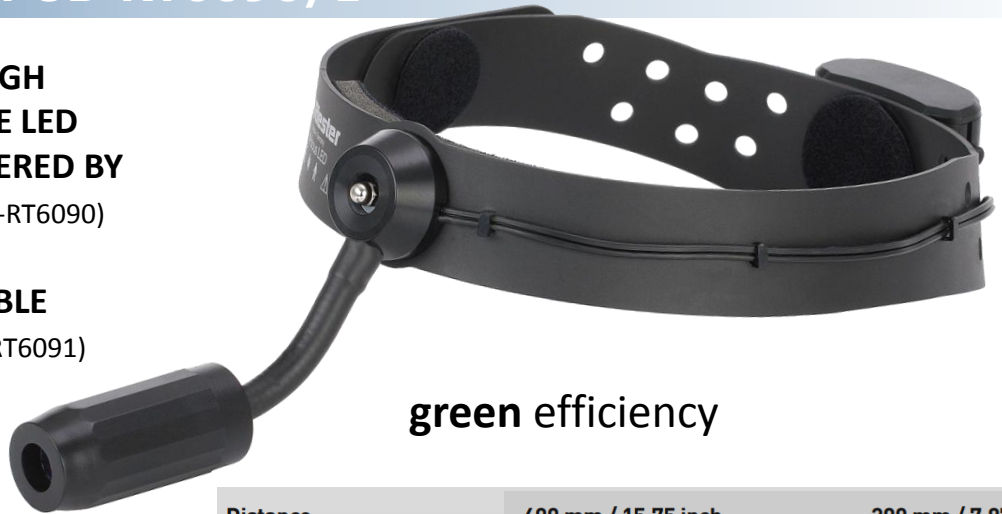

green efficiency

Distance	400 mm / 15.75 inch	200 mm / 7.87 inch
Lux max	10.000	23.000
Circle of light max. Ø	200 mm / 7.87 inch	120 mm / 4.73 inch
Circle of light min. Ø	40 mm / 1.57 inch	30 mm / 1.18 inch

DESCRIPTION:

- **One hand operation**
- **Wireless and balanced headband with built in battery compartment**
- **More efficient diagnosis due to white high performance LED (140 lumen)**
- **Service life of LED up to 50.000 hours**
- **Significantly less energy consumption and heat generation**
- Headband version, with built-in battery compartment and integral charge jack, can be optionally fitted with 2 lithium batteries size CR 123 A, 4 batteries size AAA (alkaline) or 4 rechargeable batteries size AAA (NiMH).
- Headband cleaning made easy by inner, removable and washable padding.
- Particularly comfortable with balanced, infinitely adjustable headband.
- On/off switch on battery compartment.
- Charge jack for plug-in charger integrated in battery compartment.
- Focusable LED lamp at a distance of 40 cm /15.75 inch approx. Ø 40 mm /1.57 inch - Ø 200 mm / 7.87 inch, infinitely adjustable.
- White LED with 6.500 Kelvin.
- Battery operation: approx. 4 hrs. with lithium batteries, approx. 90 min. with alkaline batteries Type AAA.
- Rechargeable battery operation: approx. 90 min. with fully charged battery.
- Weight: 210 g / 0.46 lb without batteries, 242 g / 0.54 lb with CR 123 A lithium batteries, 259 g / 0.57 lb with rechargeable AAA batteries.



TEL: 0 100 350 238

FAX : 086 477 9020

E-MAIL: info@ububelemedical.co.za